

**ACP2520-A - Oct. 28, 2015**

Item # ACP2520-A was discontinued on Oct. 28, 2015. For informational purposes, this is a copy of the website content at that time and is valid only for the stated product.

- ▶ UV, Visible, and IR Models Available
- ▶ Optimized Heat Management Results in Stable Output
- ▶ Internal SM1 (1.035"-40) Threading
- ▶ Collimation Adapters Available Separately




**M405LP1**  
 405 nm Mounted LED  
 1500 mW Output Power



**M505L3**  
 505 nm Mounted LED  
 400 mW Output Power

M625L3 with a Collimator  
 Used as a Light Source for a  
 Microscope

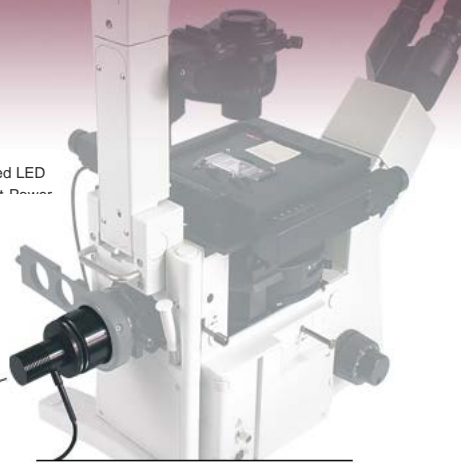

[Hide Overview](#)

OVERVIEW

**Mounted LED Features**

- Nominal Wavelengths Ranging from 265 nm to 1550 nm
- Broadband, Warm White (3000 K), and Cold White (6500 K) LEDs Also Available
- Integrated EEPROM Stores LED Operating Parameters
- Thermal Properties Optimized for Stable Output Power
- Internal SM1 (1.035"-40) Threading (6 mm Deep) for Attaching Collimation Adapters or Ø1" Lens Tubes
- Collimation Adapters Available
  - Ø2" Optics in Housings for Selected Leica, Nikon, Olympus, and Zeiss Microscopes
  - Ø1" Optics in Housings with Internally SM2-Threaded Outputs
- 4-Pin Female Mating Connector for Custom Power Supplies can be Purchased Separately
- Versions with the Collimation Adapter Included can be Found Here

Each uncollimated, mounted LED consists of a single LED that has been mounted to the end of a heat sink. Most of the LEDs on this page have a heat sink housing that has the same external diameter (1.20") as an SM1 lens tube. LEDs with a higher power output (>1000 mW) are mounted to a larger heat sink with a Ø57.0 mm plastic housing for increased heat dissipation and thermal stability. All the heat sinks are equipped with 6 mm long internal SM1 (1.035"-40) threads for easy integration with other Thorlabs components. The larger heat sinks are also equipped with four 4-40 tapped holes for compatibility with 30 mm cage systems.

The integrated EEPROM chip in each LED stores information about the LED (e.g., current limit, wavelength, and forward voltage) and can be read by Thorlabs' DC2200 and DC4100 LED Controllers. For more information about LED drivers, including the basic LEDD1B driver, see the [LED Drivers](#) tab.

**Optimized Thermal Management**

| Item #               | Color<br>(Click for Spectrum) <sup>a</sup> | Nominal Wavelength <sup>a,b</sup> | Minimum LED Power Output <sup>a</sup> |
|----------------------|--|-----------------------------------|---------------------------------------|
| M265L3 <sup>c</sup>  | Deep UV                                    | 265 nm                            | 10 mW                                 |
| M280L3 <sup>c</sup>  | Deep UV                                    | 280 nm                            | 25 mW                                 |
| M310L3 <sup>c</sup>  | Deep UV                                    | 310 nm                            | 25 mW                                 |
| M340L3 <sup>c</sup>  | UV   | 340 nm                            | 10 mW                                 |
| M365L2 <sup>c</sup>  | UV   | 365 nm                            | 190 mW                                |
| M365LP1 <sup>c</sup> | UV   | 365 nm                            | 1150 mW                               |
| M375L3 <sup>c</sup>  | UV   | 375 nm                            | 387 mW                                |
| M385L2 <sup>c</sup>  | UV   | 385 nm                            | 270 mW                                |
| M385LP1 <sup>c</sup> | UV   | 385 nm                            | 1650 mW                               |
| M395L4 <sup>c</sup>  | UV   | 395 nm                            | 400 mW                                |
| M405L2 <sup>c</sup>  | UV   | 405 nm                            | 410 mW                                |
| M405LP1 <sup>c</sup> | UV   | 405 nm                            | 1500 mW                               |
| M420L3 <sup>c</sup>  | Violet                                     | 420 nm                            | 750 mW                                |
| M450LP1              | Royal Blue                                 | 450 nm                            | 1850 mW                               |
| M455L3               | Royal Blue                                 | 455 nm                            | 900 mW                                |
| M470L3               | Blue                                       | 470 nm                            | 650 mW                                |
| M490L3               | Blue                                       | 490 nm                            | 200 mW                                |
| M505L3               | Cyan                                       | 505 nm                            | 400 mW                                |
| M530L3               | Green                                      | 530 nm                            | 350 mW                                |
| M565L3 <sup>d</sup>  | Lime                                       | 565 nm                            | 880 mW                                |
| M590L3               | Amber                                      | 590 nm                            | 160 mW                                |
| M595L3 <sup>d</sup>  | Amber                                      | 595 nm                            | 445 mW                                |
| M617L3               | Orange                                     | 617 nm                            | 600 mW                                |
| M625L3               | Red  | 625 nm                            | 700 mW                                |
| M660L4               | Deep Red                                   | 660 nm                            | 940 mW                                |
| M730L4               | Far Red                                    | 730 nm                            | 515 mW                                |
| M780L3               | IR   | 780 nm                            | 200 mW                                |
| M810L3               | IR   | 810 nm                            | 325 mW                                |
| M850L3               | IR   | 850 nm                            | 900 mW                                |
| M880L3               | IR   | 880 nm                            | 300 mW                                |
| M940L3               | IR   | 940 nm                            | 800 mW                                |
| M970L3               | IR   | 970 nm                            | 35 mW                                 |
| M1050L2              | IR   | 1050 nm                           | 50 mW                                 |

These mounted LEDs possess good thermal stability properties, eliminating the issue of degradation of optical output power due to increased LED temperature. For more details, please see the *Stability* tab.

#### Broadband LED Option

The MBB1L3 mounted LED has been designed to have relatively flat spectral emission over a wide wavelength range. Its FWHM bandwidth ranges from 500 nm to 780 nm, while the 10 dB bandwidth ranges between 470 nm and 850 nm. For more information on the spectrum of this broadband source, please see the table to the right.

#### Collimation & Microscope Adapters

Collimation adapters are available that contain an AR-coated aspheric lens for LEDs with wavelengths from 350 nm to 1050 nm. The COP adapters are designed to mate to the epillumination ports on Leica DMI, Nikon Eclipse, Olympus IX/BX, and Zeiss Axioskop microscopes; see below for more details. Additionally, Thorlabs offers mounted LEDs with microscope adapters pre-attached.

Alternatively, the ACP2520 collimation adapters

have an internally M34 x 0.5 threaded housing and use a Ø25 mm optic in a translating carriage to provide collimation adjustment. Each collimation adapter includes an external M34 x 0.5 to internal SM2 (2.035"-40) thread adapter so that the LEDs can be easily integrated with Thorlabs' SM2-threaded components and microscope port adapters. This collimation adapter is also offered without an included optic (Item # ACP), allowing a user-supplied Ø1" or Ø25 mm optic to be integrated with the LEDs. See the *Collimation* tab for more information.

#### Multi-LED Source

A customizable multi-LED source may be constructed using our mounted LEDs and other Thorlabs items. This source may be configured for integration with Thorlabs' versatile SM1 Lens Tube Systems, 30 mm Cage Systems, and the microscope adapters sold below. Please see the *Multi-LED Source* tab for a detailed item list and instructions.

Thorlabs also offers integrated, user-configurable 4-Wavelength High-Power LED Sources.

#### Driver Options

Thorlabs offers four LED drivers: LEDD1B, DC2200, DC4100, and DC4104 (the latter two require the DC4100-HUB). See the *LED Drivers* tab for compatibility and driver features. The LEDD1B is capable of providing LED modulation frequencies up to 5 kHz, while DC4100 and DC4104 can modulate the LED at a rate up to 100 kHz. The DC2200 can provide modulation at up to 250 kHz if driven by an external source. In addition, the DC2200, DC4100, and DC4104 drivers are capable of reading the current limit from the EEPROM chip of the connected LED and automatically adjusting the max current setting to protect the LED.

|                           |            |                           |        |
|---------------------------|------------|---------------------------|--------|
| <b>M1200L3</b>            | IR         | 1200 nm                   | 30 mW  |
| <b>M1300L3</b>            | IR         | 1300 nm                   | 25 mW  |
| <b>M1450L3</b>            | IR         | 1450 nm                   | 31 mW  |
| <b>M1550L3</b>            | IR         | 1550 nm                   | 31 mW  |
| <b>MBB1L3<sup>e</sup></b> | Broadband  | 470 - 850 nm <sup>f</sup> | 70 mW  |
| <b>MWWHL3<sup>d</sup></b> | Warm White | 3000 K <sup>g</sup>       | 500 mW |
| <b>MCWHL5<sup>d</sup></b> | Cold White | 6500 K <sup>g</sup>       | 800 mW |

- Due to variations in the manufacturing process and operating parameters such as temperature and current, the actual spectral output of any given LED will vary. Output plots and nominal wavelength specs are only intended to be used as a guideline.
- For LEDs in the visible spectrum, the nominal wavelength indicates the wavelength at which the LED appears brightest to the human eye. For UV and IR LEDs, the nominal wavelength corresponds to the peak wavelength. The nominal wavelength for visible LEDs may not correspond to the peak wavelength as measured by a spectrograph.
- Our 265 nm to 420 nm LEDs radiate intense UV light during operation. Precautions must be taken to prevent looking directly at the UV light and UV light protective glasses must be worn to avoid eye damage. Exposure of the skin and other body parts to the UV light should be avoided.
- The M565L3, M595L3, MWWHL3, and MCWHL5 are phosphor-converted LEDs and may not turn off completely when modulated above 10 kHz at duty cycles below 50%.
- The MBB1L3 LED may not turn off completely when modulated at frequencies above 1 kHz with a duty cycle of 50%, as the broadband emission is produced by optically stimulating emission from phosphor. For modulation at frequencies above 1 kHz, the duty cycle may be reduced. For example, 10 kHz modulation is attainable with a duty cycle of 5%.
- 10 dB bandwidth.
- Correlated color temperature.

[Hide Specs](#)

## S P E C S

Green shaded rows indicate LEDs mounted to the larger Ø57.0 mm housing.







| Item #                     | Color (Click for Spectrum) <sup>a</sup> | Nominal Wavelength <sup>a,b</sup> | Minimum LED Power Output <sup>a</sup> | Typical LED Power Output <sup>a</sup> | Maximum Current (CW) | Forward Voltage | Bandwidth (FWHM) | Irradiance (Typical) <sup>c</sup> | Electrical Power | Viewing Angle (Full Angle at Half Max) | Emitter Size             | Typical Lifetime <sup>d</sup> |
|----------------------------|---|-----------------------------------|---------------------------------------|---------------------------------------|----------------------|-----------------|------------------|-----------------------------------|------------------|--|--------------------------|-------------------------------|
| <b>M265L3<sup>e</sup></b>  | Deep UV                                 | 265 nm                            | 10 mW                                 | 12 mW                                 | 350 mA               | 6.8 V           | 11 nm            | -                                 | 2.380 W          | 130°                                   | 1 mm x 1 mm              | >1 000 h                      |
| <b>M280L3<sup>e</sup></b>  | Deep UV                                 | 280 nm                            | 25 mW                                 | 30 mW                                 | 350 mA               | 5.9 V           | 12 nm            | 3.9 µW/mm <sup>2</sup>            | 2.065 W          | 140°                                   | 1 mm x 1 mm              | >500 h                        |
| <b>M310L3<sup>e</sup></b>  | Deep UV                                 | 310 nm                            | 25 mW                                 | 30 mW                                 | 350 mA               | 5.9 V           | 10 nm            | 5.8 µW/mm <sup>2</sup>            | 2.065 W          | 140°                                   | 1 mm x 1 mm              | >500 h                        |
| <b>M340L3<sup>e</sup></b>  | UV                                      | 340 nm                            | 10 mW                                 | 12 mW                                 | 80 mA                | 8.1 V           | 8 nm             | 0.1 µW/mm <sup>2</sup>            | 0.648 W          | 176°                                   | 2 mm x 2 mm <sup>f</sup> | >3 000 h                      |
| <b>M365L2<sup>e</sup></b>  | UV                                      | 365 nm                            | 190 mW                                | 360 mW                                | 700 mA               | 4.4 V           | 7.5 nm           | 8.9 µW/mm <sup>2</sup>            | 3.080 W          | 120°                                   | 1 mm x 1 mm              | >10 000 h                     |
| <b>M365LP1<sup>e</sup></b> | UV                                      | 365 nm                            | 1150 mW                               | 1400 mW                               | 1400 mA              | 3.75 V          | 9 nm             | 17.6 µW/mm <sup>2</sup>           | 5.250 W          | 120°                                   | 1.4 mm x 1.4 mm          | >10 000 h                     |
| <b>M375L3<sup>e</sup></b>  | UV                                      | 375 nm                            | 387 mW                                | 470 mW                                | 700 mA               | 3.8 V           | 9 nm             | 14.1 µW/mm <sup>2</sup>           | 2.660 W          | 110°                                   | 1 mm x 1 mm              | >10 000 h                     |
| <b>M385L2<sup>e</sup></b>  | UV                                      | 385 nm                            | 270 mW                                | 430 mW                                | 700 mA               | 4.3 V           | 10 nm            | 11.8 µW/mm <sup>2</sup>           | 3.010 W          | 120°                                   | 1 mm x 1 mm              | >10 000 h                     |
| <b>M385LP1<sup>e</sup></b> | UV                                      | 385 nm                            | 1650 mW                               | 1830 mW                               | 1400 mA              | 3.65 V          | 12 nm            | 23.3 µW/mm <sup>2</sup>           | 5.110 W          | 120°                                   | 1.4 mm x 1.4 mm          | >10 000 h                     |
| <b>M395L4<sup>e</sup></b>  | UV                                      | 395 nm                            | 400 mW                                | 535 mW                                | 500 mA               | 4.5 V           | 16 nm            | 6.7 µW/mm <sup>2</sup>            | 2.250 W          | 126°                                   | 1 mm x 1 mm              | >10 000 h                     |
| <b>M405L2<sup>e</sup></b>  | UV                                      | 405 nm                            | 410 mW                                | 760 mW                                | 1000 mA              | 3.8 V           | 13 nm            | 37.1 µW/mm <sup>2</sup>           | 3.800 W          | 85°                                    | 1 mm x 1 mm              | 100 000 h                     |
| <b>M405LP1<sup>e</sup></b> | UV                                      | 405 nm                            | 1500 mW                               | 1700 mW                               | 1400 mA              | 3.45 V          | 12 nm            | 24.6 µW/mm <sup>2</sup>           | 4.830 W          | 120°                                   | 1.4 mm x 1.4 mm          | >10 000 h                     |

|                     |            |                           |         |         |         |        |        |                                 |         |      |                 |           |
|---------------------|------------|---------------------------|---------|---------|---------|--------|--------|---------------------------------|---------|------|-----------------|-----------|
| M420L3 <sup>e</sup> | Violet     | 420 nm                    | 750 mW  | 820 mW  | 1000 mA | 3.5 V  | 15 nm  | 13.1 $\mu\text{W}/\text{mm}^2$  | 3.500 W | 125° | 1 mm x 1 mm     | >10 000 h |
| M450LP1             | Royal Blue | 450 nm                    | 1850 mW | 2100 mW | 2000 mA | 3.5 V  | 18 nm  | 35.6 $\mu\text{W}/\text{mm}^2$  | 7.000 W | 120° | 1.5 mm x 1.5 mm | 1 000 h   |
| M455L3              | Royal Blue | 455 nm                    | 900 mW  | 1020 mW | 1000 mA | 3.2 V  | 18 nm  | 31.2 $\mu\text{W}/\text{mm}^2$  | 3.200 W | 80°  | 1 mm x 1 mm     | 100 000 h |
| M470L3              | Blue       | 470 nm                    | 650 mW  | 710 mW  | 1000 mA | 3.2 V  | 25 nm  | 21.9 $\mu\text{W}/\text{mm}^2$  | 3.200 W | 80°  | 1 mm x 1 mm     | 100 000 h |
| M490L3              | Blue       | 490 nm                    | 200 mW  | 250 mW  | 350 mA  | 3.5 V  | 23 nm  | 15.7 $\mu\text{W}/\text{mm}^2$  | 1.225 W | 22°  | 1 mm x 1 mm     | >10 000 h |
| M505L3              | Cyan       | 505 nm                    | 400 mW  | 440 mW  | 1000 mA | 3.3 V  | 30 nm  | 11.1 $\mu\text{W}/\text{mm}^2$  | 3.300 W | 80°  | 1 mm x 1 mm     | 100 000 h |
| M530L3              | Green      | 530 nm                    | 350 mW  | 370 mW  | 1000 mA | 3.2 V  | 33 nm  | 9.5 $\mu\text{W}/\text{mm}^2$   | 3.200 W | 80°  | 1 mm x 1 mm     | 100 000 h |
| M565L3 <sup>g</sup> | Lime       | 565 nm                    | 880 mW  | 979 mW  | 1000 mA | 3.1 V  | 104 nm | 11.7 $\mu\text{W}/\text{mm}^2$  | 3.100 W | 125° | 1 mm x 1 mm     | 50 000 h  |
| M590L3              | Amber      | 590 nm                    | 160 mW  | 170 mW  | 1000 mA | 2.2 V  | 18 nm  | 5.3 $\mu\text{W}/\text{mm}^2$   | 2.200 W | 80°  | 1 mm x 1 mm     | 100 000 h |
| M595L3 <sup>g</sup> | Amber      | 595 nm                    | 445 mW  | 502 mW  | 700 mA  | 3.05 V | 80 nm  | 6.9 $\mu\text{W}/\text{mm}^2$   | 2.135 W | 120° | 1 mm x 1 mm     | 50 000 h  |
| M617L3              | Orange     | 617 nm                    | 600 mW  | 650 mW  | 1000 mA | 2.2 V  | 18 nm  | 15.7 $\mu\text{W}/\text{mm}^2$  | 2.200 W | 80°  | 1 mm x 1 mm     | 100 000 h |
| M625L3              | Red        | 625 nm                    | 700 mW  | 770 mW  | 1000 mA | 2.2 V  | 18 nm  | 18.0 $\mu\text{W}/\text{mm}^2$  | 2.200 W | 80°  | 1 mm x 1 mm     | 100 000 h |
| M660L4              | Deep Red   | 660 nm                    | 940 mW  | 1050 mW | 1200 mA | 2.6 V  | 20 nm  | 20.88 $\mu\text{W}/\text{mm}^2$ | 3.120 W | 120° | 1.5 mm x 1.5 mm | >10 000 h |
| M730L4              | Far Red    | 730 nm                    | 515 mW  | 595 mW  | 1000 mA | 2.3 V  | 37 nm  | 13.2 $\mu\text{W}/\text{mm}^2$  | 2.300 W | 160° | 1 mm x 1 mm     | >10 000 h |
| M780L3              | IR         | 780 nm                    | 200 mW  | 300 mW  | 800 mA  | 2.0 V  | 28 nm  | 47.3 $\mu\text{W}/\text{mm}^2$  | 1.600 W | 20°  | 1 mm x 1 mm     | >10 000 h |
| M810L3              | IR         | 810 nm                    | 325 mW  | 375 mW  | 500 mA  | 3.6 V  | 25 nm  | 61.8 $\mu\text{W}/\text{mm}^2$  | 1.800 W | 40°  | 1 mm x 1 mm     | >10 000 h |
| M850L3              | IR         | 850 nm                    | 900 mW  | 1100 mW | 1000 mA | 2.9 V  | 30 nm  | 22.9 $\mu\text{W}/\text{mm}^2$  | 2.900 W | 90°  | 1 mm x 1 mm     | 100 000 h |
| M880L3              | IR         | 880 nm                    | 300 mW  | 350 mW  | 1000 mA | 1.7 V  | 50 nm  | 5.6 $\mu\text{W}/\text{mm}^2$   | 1.700 W | 128° | 1 mm x 1 mm     | >10 000 h |
| M940L3              | IR         | 940 nm                    | 800 mW  | 1000 mW | 1000 mA | 2.75 V | 37 nm  | 19.1 $\mu\text{W}/\text{mm}^2$  | 2.750 W | 90°  | 1 mm x 1 mm     | 100 000 h |
| M970L3              | IR         | 970 nm                    | 35 mW   | 50 mW   | 600 mA  | 1.4 V  | 50 nm  | 0.7 $\mu\text{W}/\text{mm}^2$   | 0.840 W | 124° | 1 mm x 1 mm     | >10 000 h |
| M1050L2             | IR         | 1050 nm                   | 50 mW   | 70 mW   | 700 mA  | 1.5 V  | 60 nm  | 1.9 $\mu\text{W}/\text{mm}^2$   | 1.050 W | 120° | 1 mm x 1 mm     | >10 000 h |
| M1200L3             | IR         | 1200 nm                   | 30 mW   | 35 mW   | 700 mA  | 1.4 V  | 80 nm  | 0.7 $\mu\text{W}/\text{mm}^2$   | 0.980 W | 134° | 1 mm x 1 mm     | >10 000 h |
| M1300L3             | IR         | 1300 nm                   | 25 mW   | 30 mW   | 500 mA  | 1.4 V  | 80 nm  | 0.6 $\mu\text{W}/\text{mm}^2$   | 0.700 W | 134° | 1 mm x 1 mm     | >10 000 h |
| M1450L3             | IR         | 1450 nm                   | 31 mW   | 36 mW   | 700 mA  | 1.15 V | 80 nm  | 0.4 $\mu\text{W}/\text{mm}^2$   | 0.805 W | 136° | 1 mm x 1 mm     | >10 000 h |
| M1550L3             | IR         | 1550 nm                   | 31 mW   | 36 mW   | 700 mA  | 1.5 V  | 102 nm | 0.5 $\mu\text{W}/\text{mm}^2$   | 1.050 W | 136° | 1 mm x 1 mm     | >10 000 h |
| MBB1L3 <sup>h</sup> | Broadband  | 470 - 850 nm <sup>i</sup> | 70 mW   | 80 mW   | 500 mA  | 3.6 V  | 280 nm | 12.5 $\mu\text{W}/\text{mm}^2$  | 1.800 W | 120° | 1 mm x 1 mm     | 10 000 h  |
| MWWHL3 <sup>9</sup> | Warm White | 3000 K <sup>j</sup>       | 500 mW  | 550 mW  | 1000 mA | 3.1 V  | N/A    | 10.7 $\mu\text{W}/\text{mm}^2$  | 3.100 W | 120° | 1 mm x 1 mm     | >50 000 h |
| MCWHL5 <sup>9</sup> | Cold White | 6500 K <sup>j</sup>       | 800 mW  | 840 mW  | 1000 mA | 3.2 V  | N/A    | 24.8 $\mu\text{W}/\text{mm}^2$  | 3.200 W | 80°  | 1 mm x 1 mm     | 100 000 h |

- Due to variations in the manufacturing process and operating parameters such as temperature and current, the actual spectral output of any given LED will vary. Output plots and nominal wavelength specs are only intended to be used as a guideline.
- For LEDs in the visible spectrum, the nominal wavelength indicates the wavelength at which the LED appears brightest to the human eye. For UV and IR LEDs, the nominal wavelength corresponds to the peak wavelength. The nominal wavelength for visible LEDs may not correspond to the peak wavelength as measured by a spectrograph.
- Irradiance is measured at a distance of 200 mm from the LED.
- Thorlabs defines the lifetime of our LEDs as  $B_{50}/L_{50}$ , meaning that 50% of the LEDs with a given item # will fall below 50% of the initial optical power at the end of the specified lifetime. Please see the *Stability* tab for more details.
- Our 265 nm to 420 nm LEDs radiate intense UV light during operation. Precautions must be taken to prevent looking directly at the UV light and UV light protective glasses must be worn to avoid eye damage. Exposure of the skin and other body parts to the UV light should be avoided.
- The M340L3 LED is comprised of four 1 mm x 1 mm emitters that have a total area of 2 mm x 2 mm.
- The M565L3, M595L3, MWWHL3, and MCWHL5 are phosphor-converted LEDs and may not turn off completely when modulated above 10 kHz at duty cycles below 50%.
- The MBB1L3 LED may not turn off completely when modulated at frequencies above 1 kHz with a duty cycle of 50%, as the broadband emission is produced by optically stimulating emission from phosphor. For modulation at frequencies above 1 kHz, the duty cycle may be reduced. For example, 10 kHz modulation is attainable with a duty cycle of 5%.
- 10 dB Bandwidth.
- Correlated Color Temperature.

### Relative Power

The actual spectral output and total output power of any given LED will vary due to variations in the manufacturing process and operating parameters, such as temperature and current. Both a typical and minimum output power are specified to help you select an LED that suits your needs. Each mounted high-power LED will provide at least the minimum specified output power at the maximum current. In order to provide a point of comparison for the relative powers of LEDs with different nominal wavelengths, the spectra in the plots below have been scaled to the minimum output power for each LED. This data is representative, not absolute. An excel file with normalized and scaled spectra for all of the mounted high-power LEDs can be downloaded [here](#).

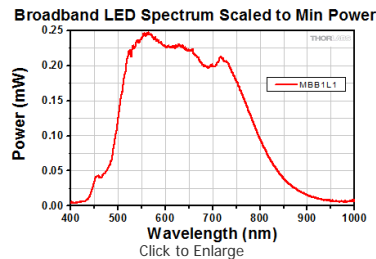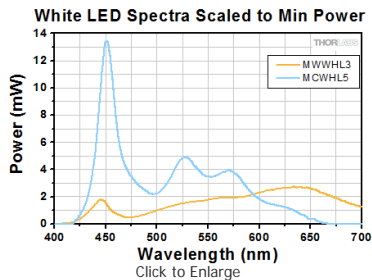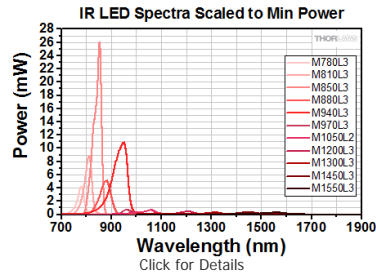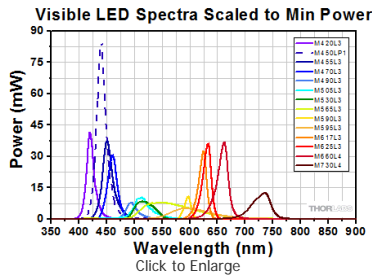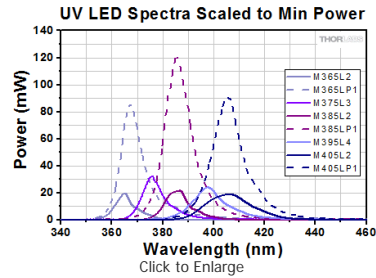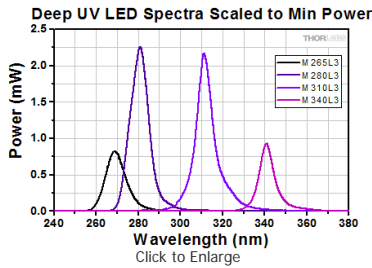


[Hide Stability](#)

### STABILITY


#### LED Lifetime and Long-Term Power Stability

One characteristic of LEDs is that they naturally exhibit power degradation with time. Often this power degradation is slow, but there are also instances where large, rapid drops in power, or even complete LED failure, occur. LED lifetimes are defined as the time it takes a specified percentage of a type of LED to fall below some power level. The parameters for the lifetime measurement can be written using the notation  $B_{XX}/L_{YY}$ , where XX is the percentage of that type of LED that will provide less than YY percent of the specified output power after the lifetime has elapsed. Thorlabs defines the lifetime of our LEDs as  $B_{50}/L_{50}$ , meaning that 50% of the LEDs with a given Item # will fall below 50% of the initial optical power at the end of the specified lifetime. For example, if a batch of 100 LEDs is rated for 150 mW of output power, 50 of these LEDs can be expected to produce an output power of  $\leq 75$  mW after the specified LED lifetime has elapsed.


The sample plots below show example data from long-term stability testing for our UV LEDs over a 45 day period; the 280 nm and 310 nm LEDs have a typical lifetime of >500 hours (~20 days), while the M340L3 has a lifetime of >3,000 hours (~125 days). The small power drop experienced by each LED after it is turned on is typical behavior during the first few minutes of operation. It corresponds to the period of time required for the LED to warm up to the point where it is thermally stable. Please note that each graph represents the performance of a single LED; performance of individual LEDs will vary within the stated specifications.

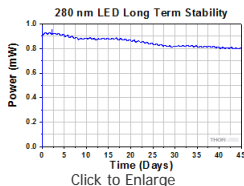

The M280L3 LED has a typical lifetime of >500 hours. In this case, the unit under test had dropped to ~80% of the initial output power after 45 days.

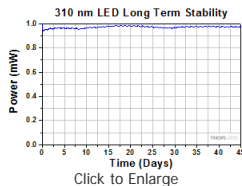


The M310L3 LED has a typical lifetime of >500 hours. In this case, the unit under test continued to provide more than 95% of its initial power after 45 days.

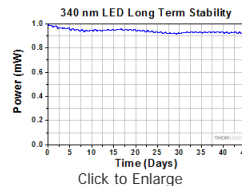

The M340L3 LED has a typical lifetime of >3,000 hours. In this case, the unit under test continued to provide more than 90% of its initial power after 45 days.

#### Optimized Thermal Management

The thermal dissipation performance of these mounted LEDs has been optimized for stable power output. The heat sink is directly mounted to the LED mount so as to provide optimal thermal contact. By doing so, the degradation of optical output power that can be attributed to increased LED junction temperature is minimized (see the graph to the right).

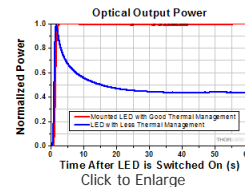

PIN DIAGRAM

Pin Connection - Male

The diagram to the right shows the male connector of the mounted LED assembly. It is a standard M8 x 1 sensor circular connector. Pins 1 and 2 are the connection to the LED. Pin 3 and 4 are used for the internal EEPROM in these LEDs. If using an LED driver that was not purchased from Thorlabs, be careful that the appropriate connections are made to Pin 1 and Pin 2 and that you do not attempt to drive the LED through the EEPROM pins.



| Pin | Specification | Color |
|-----|---------------|-------|
| 1   | LED Anode     | Brown |
| 2   | LED Cathode   | White |
| 3   | EEPROM GND    | Black |
| 4   | EEPROM IO     | Blue  |

LED DRIVERS

| Compatible Drivers                                     | LEDD1B <sup>a</sup>   | DC2200 <sup>b</sup>   | DC4100 <sup>b,c,d</sup>                                    | DC4104 <sup>b,c,d</sup>                                     |
|--|---|---|--|---|
| Click Photos to Enlarge                                |   |   |  |   |
| Max LED Driver Current Output                          | 1.2 A   | LED1 Terminal: 10.0 A<br>LED2 Terminal: 2.0 A <sup>e</sup>                                      | 1.0 A per Channel  | 1.0 A per Channel   |
| Max LED Driver Forward Voltage                         | 12 V  | 50 V  | 5 V  | 5 V   |
| Max Modulation Frequency Using External Input          | 5 kHz   | 250 kHz <sup>f,g</sup>  | 100 kHz <sup>g</sup><br>(Simultaneous Across all Channels) | 100 kHz <sup>g</sup><br>(Independently Controlled Channels) |
| External Control Interface(s)                          | Analog (BNC)  | USB 2.0 and Analog (BNC)  | USB 2.0 and Analog (BNC)                                   | USB 2.0 and Analog (8-Pin)                                  |
| Main Driver Features                                   | Very Compact Footprint<br>60 mm x 73 mm x 104 mm<br>(W x H x D) | Touchscreen Interface with Internal and External Options for Pulsed and Modulated LED Operation | 4 Channels <sup>c</sup>                                    | 4 Channels <sup>c</sup>                                     |
| EEPROM Compatible: Reads Out LED Data for LED Settings | -   | ✓   | ✓  | ✓   |
| LCD Display  | -   | ✓   | ✓  | ✓   |

- The LEDD1B should not be used to drive the M340L3, as the current limit can only be set to a minimum of 200 mA (compared to the M340L3's max drive current of 80 mA).
- Automatically limits to LED's max current via EEPROM readout.
- The DC4100 and DC4104 can power and control up to four LEDs simultaneously when used with the DC4100-HUB. The LEDs on this page all require the DC4100-HUB when used with the DC4100 or DC4104.
- These LED drivers have a maximum forward voltage rating of 5 V and can provide a maximum current of 1000 mA. As a result, they cannot be used to drive the M280L3, M310L3, or M340L3 LEDs, which have forward voltage ratings of >5 V. They can be used to drive the M660L3 LED, but will not be able to provide the LED's maximum current of 1200 mA.
- The mounted LEDs sold below are compatible with the LED2 Terminal.
- Small Signal Bandwidth: Modulation not exceeding 20% of full scale current. The driver accepts other waveforms, but the maximum frequency will be reduced.
- Several of these LEDs produce light by stimulating emission from phosphor, which limits their modulation frequencies. The M565L3, M595L3, MWWHL3, and MCWHL5 LEDs may not turn off completely when modulated above 10 kHz at duty cycles below 50%. The MBB1L3 LED may not turn off completely when modulated at frequencies above 1 kHz with a duty cycle of 50%. When the MBB1L3 is modulated at frequencies above 1 kHz, the duty cycle may be reduced; for example, 10 kHz modulation is attainable with a duty cycle of 5%.

Note: The DC3100 drivers sold with our Modulated LEDs for FLIM Microscopy kits are not compatible with the LEDs sold on this page.

COLLIMATION

Collimating the LED

Thorlabs' extensive catalog of mechanical and optical components provides a variety of configurations that can be used to collimate our mounted LEDs. Some of the applications of the collimated LEDs include custom imaging systems, microscope illuminators, or projectors. Our COP adapters, available below, feature microscope-compatible outputs and Ø2" aspheric condenser lenses. The ACP, ACP2520-A, and ACP2520-B, also available below, are designed to provide collimation with focus adjustment via a Ø1" or Ø25 mm optic in a translating carriage; the mechanical housings have an M34 x 0.5 internal output thread and include an internally SM2-threaded adapter for integration with standard Thorlabs' components. If your setup requires a collimation package with the smallest possible profile, the LEDs can also be integrated with Ø1" collimating optics and SM1-threaded lens tubes. When exchanging the lens in your collimation adapter, please be careful to follow proper optics handling procedures (Optic Handling Tutorial).

ACP Collimation Adapter

Thorlabs' ACP collimation adapters accept Ø1" or Ø25 mm collimation optics. The ACP2520-A and ACP2520-B include an aspheric condenser lenses coated for 350 - 700 nm (Item # ACL2520U-A) and 650 - 1050 nm (Item # ACL2520U-B), respectively. The ACP is offered without a pre-installed optics so that user-supplied components can be integrated with the LEDs. Several suggestions are presented in the table to the right.

Installing a new lens in the ACP is a simple procedure:

| Item #  | Qty | Description  |
|---|-----|--|
| SPW602 <sup>a</sup>   | 1   | Spanner Wrench for SM1-Threaded Retaining Rings    |
| ACL2520U-DG6-A <sup>b</sup> ,<br>ACL2520U-DG6-B <sup>b</sup> ,<br>ACL2520U <sup>b</sup> ,<br>or ACL2520U-DG6 <sup>b</sup> | 1   | Aspheric Condenser Lens (with or without Diffuser) |

- The ACP includes a retaining ring that is thicker than our standard



Click to Enlarge  
An ACP2520-A Collimation  
Adapter Installed on a  
Mounted LED

1. Turn the adjustment knob to move the optic mounting carriage to the M34 x 0.5 threaded end of the housing.
2. Use the SPW602 spanner wrench to remove the SM1-threaded retaining ring from the housing.
3. Place a Ø1" or Ø25 mm optic of your choice into the mounting carriage with the curved surface facing the output. For customers concerned with the homogeneity of the beam, the AR-coated aspheric condenser lens with diffuser (ACL2520U-DG6-A or ACL2520U-DG6-B) is a good option.
4. Use the SPW602 spanner wrench to screw the SM1-threaded retaining ring into the mounting carriage, securing the optic in place.
5. Screw the externally SM1-threaded end of the collimation adapter onto an LED of your choice as shown in the picture to the left.

SM1RR SM1-threaded retaining ring so that the SPW602 can be used without scratching the highly curved surface of an aspheric condenser lens.

- -A and -B refer to the type of AR coating on the lens. Thorlabs' LEDs with a nominal wavelength between 365 nm and 660 nm would require the -A coating, while the LEDs with a nominal wavelength between 730 nm and 1050 nm would require the -B Coating. IR LEDs that emit past 1050 nm can be collimated using an uncoated condenser lens, such as the ACL2520U. Deep UV LEDs with a wavelength  $\leq 340$  nm require a lens fabricated from UV Fused Silica, since many standard varieties of glass do not transmit below 350 nm.

### SM1-Threaded Collimation Assembly

For cases where a smaller profile than the ACP collimation adapter is required, a simple LED collimation assembly can be built from the components listed in the table to the lower right.

1. First, install the optic in the adjustable lens tube, which allows one to control the working distance of the lens while collimating the LED. The SM1-threaded (1.035"-40) SM1V05 comes with a locking nut and a retaining ring. For customers concerned with the homogeneity of the beam, the AR-coated aspheric condenser lens with diffuser (ACL2520U-DG6-A or ACL2520U-DG6-B) is a good option. By the end of this step, the lens will rest on top of one retaining ring (SM1RR) and be secured in place by another retaining ring placed on top of it.

| Item #   | Qty | Description  |
|--|-----|--|
| SM1RR  | 2   | Ø1" Retaining Ring<br>(One Each Included with SM1V05 & SM1L03) |
| SPW801   | 1   | Adjustable Spanner Wrench <sup>a</sup>                         |
| ACL2520U-A <sup>b</sup> , ACL2520U-B <sup>b</sup> ,<br>ACL2520U-DG6-A <sup>b</sup> , ACL2520U-DG6-B <sup>b</sup> ,<br>ACL2520U <sup>b</sup> , or ACL2520U-DG6 <sup>b</sup> | 1   | Aspheric Condenser Lens<br>(with or without Diffuser)          |
| SM1V05 <sup>c</sup>  | 1   | Ø1" Rotating Adjustable Length<br>Lens Tube, 1/2" Long         |
| SM1L03   | 1   | Ø1" Lens Tube, 0.30" Long                                      |

- a. Use the spanner wrench (SPW801) to turn the included retaining ring in the adjustable length lens tube so that it is closer to the inside lip of the tube.
- b. Carefully place the lens inside the adjustable length lens tube with the curved side facing away from the male-threaded end of the tube.
- c. Secure the lens in place with another retaining ring (SM1RR)

- While these components are SM1 threaded, we recommend our adjustable spanner wrench due to the steep curvature of the aspheric condenser lens.
- -A and -B refer to the type of AR coating on the lens. Thorlabs' LEDs with a nominal wavelength between 365 nm and 660 nm would require the -A coating, while the LEDs with a nominal wavelength between 730 nm and 1050 nm would require the -B Coating. IR LEDs that emit past 1050 nm can be collimated using an uncoated condenser lens, such as the ACL2520U. Deep UV LEDs with a wavelength  $\leq 340$  nm require a lens fabricated from UV Fused Silica, since many standard varieties of glass do not transmit below 350 nm.
- The SM1V10 Adjustable Lens Tube can also be used for this application, however, the translation range of the optic cell will be reduced from 7.6 mm (thread depth of the SM1L03) to 6 mm (thread depth of the LED). If used, the SM1L03 would no longer be needed in the assembly.

2. Thread the male end of the SM1L03 lens tube into the female end of the LED and gently tighten it.
3. Partially thread the male end of the SM1V05 adjustable length lens tube assembly into the female end of the SM1L03-LED assembly.



Click to Enlarge  
Setup for Adjustable Length Lens  
Tube and Lens



Click to Enlarge  
Completed Assembly



Click to Enlarge  
Setup for Lens Tubes and LED Assembly

### Obtaining a Well-Collimated Beam

After installing the chosen mounting adapter on a mounted LED, the distance of the lens from the LED should be adjusted by following the steps below.

A well-collimated beam has minimal divergence and will not converge at any point in the beam path. Be advised that due to the nature of the output from the LED (high emitter surface area), the beam cannot be perfectly collimated. Please refer to the table below for divergence data.

1. Power on the LED and check to see if it is properly collimated. It is easiest to check that the beam is collimated by noting the changes in the beam diameter over a range of about 1" to 2 feet away; change the distance of the lens from the LED and check again. Do this until the least divergent, non-converging, homogenous beam is obtained. The beam should be somewhat circular in diameter, may have a slightly polygonal shape, and should not be a clear image of the LED itself.
2. If you see an image of the LED, this means that the lens is not close enough to the LED. Move the lens closer to the LED until the image blurs and becomes homogenous – this is the point of collimation. Note: If the lens needs to be closer to the LED when using the SM1V05 assembly, use only one retaining ring to secure the lens in the SM1V05 so that the lens will rest on the inside lip of the SM1V05 adjustable length lens tube.





Click to Enlarge  
Image of the LED



Click to Enlarge  
Uncollimated Beam



Click to Enlarge  
Collimated Beam

5. Once the proper collimation position of the lens has been found, lock the position of the lens in place.
  - a. For the ACP collimation adapters, simply tighten the locking screw using a 2 mm (5/64") hex key.
  - b. For the SM1V05 assembly described above, loosen it from the SM1L03 lens tube by about ¼ to ½ turn, rotate the external locking nut until it is flush with the edge of the SM1L03 lens tube, and gently tighten both the assembly and the locking nut by ¼ to ½ turn (there should be slight resistance; do not over tighten). This will lock the collimation position in place.

The table below provides examples of how the half viewing angle changes for select LEDs with the addition of a Ø1" aspheric condenser lens.

| Item # | Color      | Nominal Wavelength <sup>a</sup> | Optimum Lens to Emitter Distance <sup>b</sup> | Half Viewing Angle <sup>c</sup> |                              |                                 |
|--------|------------|---------------------------------|---|---------------------------------|------------------------------|---------------------------------|
|        |            |                                 |   | +1 mm Out of Focus <sup>d</sup> | at Optimum Focusing Distance | -1 mm Out of Focus <sup>d</sup> |
| M365L2 | UV         | 365 nm                          | 12.7 mm                                       | 2.79°                           | 1.32°                        | 3.11°                           |
| M385L2 | UV         | 385 nm                          | 12.8 mm                                       | 2.68°                           | 1.33°                        | 3.06°                           |
| M405L2 | UV         | 405 nm                          | 12.9 mm                                       | 2.94°                           | 1.63°                        | 3.06°                           |
| M505L3 | Cyan       | 505 nm                          | 13.2 mm                                       | 3.52°                           | 2.72°                        | 3.46°                           |
| M625L3 | Red        | 625 nm                          | 14.4 mm                                       | 3.46°                           | 2.27°                        | 3.13°                           |
| M850L3 | IR         | 850 nm                          | 13.8 mm                                       | 3.29°                           | 3.10°                        | 3.93°                           |
| M940L3 | IR         | 940 nm                          | 13.9 mm                                       | 3.42°                           | 2.46°                        | 3.70°                           |
| MCWHL5 | Cold White | 6500 K <sup>e</sup>             | 13.9 mm                                       | 3.41°                           | 2.47°                        | 3.14°                           |

- The specifications listed in the table above are nominal values specified by the LED manufacturer.
- Optimum distance between the respective mounted LED and the ACL2520U lens used to collimate the beam.
- Power loss to  $1/e^2$  (13.5%).
- ±1 mm out of focus from Optimum Distance between the respective mounted LED and the ACL2520U lens used to collimate the beam.
- Correlated Color Temperature.

The divergence data was calculated using Zemax.

[Hide Multi-LED Source](#)

## MULTI-LED SOURCE

### Creating a Custom Multi-LED Source for Microscope Illumination

Thorlabs offers the items necessary to create your own custom multi-LED light source using two or three of the mounted LEDs offered below. As configured in the following example, the light source is intended to be used with the illumination port of a microscope. However, it may be integrated with other applications using Thorlabs' versatile SM1 Lens Tube and 30 mm Cage Systems. Thorlabs also offers integrated, user-configurable 4-Wavelength LED Sources.

#### Design & Construction

First, light will be collimated by lenses mounted in lens tubes. Dichroic mirrors mounted in kinematic cage cubes then combine the output from the multiple LEDs. The mounted LEDs may be driven by LEDD1B Compact T-Cube LED Drivers (TPS001 T-Cube Power Supplies are sold separately). The LEDD1B LED Drivers allow each LED's output to be independently modulated and can provide up to 1200 mA of current. Please take care not to drive the LED sources above their max current ratings.

When designing your custom source, select mounted LEDs from below along with dichroic mirror(s) that have cutoff wavelength(s) between the LED wavelengths. The appropriate dichroic mirror(s) will reflect light from side-mounted LEDs and transmit light along the optical axis. Please note that most of these dichroic mirrors are "longpass" filters, meaning they transmit the longer wavelengths and reflect the shorter wavelengths. To superimpose light from three or more LEDs, add each in series (as shown below), starting from the back with longer wavelength LEDs when using longpass filters. Shortpass filters may also be used if the longer wavelength is reflected and the shorter wavelength is transmitted. Sample combinations of compatible dichroic mirrors and LEDs are offered in the three tables below.

It is also necessary to select an aspheric condenser lens for each source with AR coatings appropriate for the source. Before assembling the light source, collimate the light from each mounted LED as detailed in the *Collimation* tab. For mounting the aspheric lenses in the SM1V05 Lens Tubes using the included SM1RR retaining rings, we recommend the SPW801 Adjustable Spanner Wrench. A properly collimated LED source should have a resultant beam that is approximately homogenous and not highly divergent at a distance of approximately 2 feet (60 cm). An example of a well-collimated beam is shown on the *Collimation* tab.

After each LED source is collimated, thread the SM1V05 Lens Tubes at the end of each collimated LED assembly into their respective C4W Cage Cube ports using SM1T2 Lens Tube Couplers. Install each dichroic filter in an FFM1 Dichroic Filter Holder, and mount each filter holder onto a B4C Kinematic Cage Cube Platform. Each platform is then installed in the C4W Cage Cubes by partially threading the included screws into the bottom of the cube, and then inserting and rotating the B4C platform into place. Align the platform to the desired position and then firmly tighten the screws. To connect multiple cage cubes and the microscope adapter, use the remaining SM1T2 lens tube couplers along with an SM1L05 0.5" Lens Tube between adjacent cage cubes. Finally, adjust the rotation, tip, and tilt of each B4C platform to align the reflected and transmitted beams so they overlap as closely as possible.

If desired, a multi-LED source may be constructed that employs more than three LEDs. The limiting factors for the number of LEDs that can be practically used are the collimation of the light and the dichroic mirror efficiency over the specified range. Heavier multi-LED sources may be supported with our Ø1" or Ø1.5" Posts.



Click to Enlarge  
Multi-LED Source Coupled to Microscope  
Illumination Port



Click to Enlarge  
Three-LED Source Using Components Mounted LEDs and Dichroic Mirrors  
Detailed in Example Configuration 1



Click to Enlarge  
Beam Profile of Source with 3  
Mounted LEDs



Click to Enlarge  
Two-LED source. This is  
the same as Example 1,  
but with the blue LED  
removed.

| Parts List |  |                         |           |        |
|------------|--|-------------------------|-----------|--------|
| #          | Product Description                                  | Item #                  | 2 LEDs    | 3 LEDs |
|            |  |                         | Item Qty. |        |
| 1          | Microscope Illumination Port Adapter:                | Olympus IX or BX        | 1         | 1      |
|            |  | Leica DMI               |           |        |
|            |  | Zeiss Axioskop          |           |        |
|            |  | Nikon Eclipse Ti        |           |        |
| 2          | Mounted LED <sup>b</sup>                             | -                       | 2         | 3      |
| -          | T-Cube LED Driver, 1200 mA Max Drive Current         | LEDD1B <sup>c</sup>     | 2         | 3      |
| -          | 15 V Power Supply Unit for a Single T-Cube           | TPS001 <sup>c</sup>     | 2         | 3      |
| 3          | 4-Way Mounting 30 mm Cage Cube                       | C4W                     | 1         | 2      |
| 4          | Kinematic Cage Cube Platform for C4W/C6W             | B4C                     | 1         | 2      |
| 5          | 30 mm Cage-Compatible Dichroic Filter Mount          | FFM1                    | 1         | 2      |
| 6          | Dichroic Filter(s) <sup>d</sup>                      | -                       | 1         | 2      |
| 7          | Externally SM1-Threaded End Cap                      | SM1CP2                  | 1         | 2      |
| 8          | SM1 (1.035"-40) Coupler, External Threads, 0.5" Long | SM1T2                   | 3         | 5      |
| 9          | Ø1" SM1 Lens Tube, 1/2" Long External Threads        | SM1V05                  | 2         | 3      |
| -          | Aspheric Condenser Lens                              | AR-Coated 350 - 700 nm  | 2         | 3      |
|            |  | AR-Coated 650 - 1050 nm |           |        |
| 10         | SM1 Lens Tube, 0.3" Thread Depth                     | SM1L03                  | 2         | 4      |
| -          | Blank Cover Plate with Rubber O-Ring for C4W/C6W     | B1C <sup>c</sup>        | 1         | 2      |

- The SM1A23 Zeiss Axioskop Microscope Adapter is shown.
- Mounted LEDs are available below.
- Item not pictured.
- Please see the following tables for suggested compatible LED and dichroic filter combinations, or create your own by taking into account the transmission and reflection wavelength ranges of our Dichroic Filters.
- Lenses are mounted in the SM1V05 Lens Tube in front of each LED. For each lens, select an AR coating corresponding to the emission wavelength of the LED source.

| Example Configuration 1 |          | Example Configuration 2 |          | Example Configuration 3 |          |
|-------------------------|----------|-------------------------|----------|-------------------------|----------|
| Mounted LEDs            |          | Mounted LEDs            |          | Mounted LEDs            |          |
| #                       | Item #   | #                       | Item #   | #                       | Item #   |
| 2a                      | M625L3   | 2a                      | M625L3   | 2a                      | M1050L2  |
| 2b                      | M530L3   | 2b                      | M455L3   | 2b                      | MCWHL5   |
| 2c                      | M455L3   | 2c                      | M1050L2  | Dichroic Filter(s)      |          |
| Dichroic Filter(s)      |          | Dichroic Filter(s)      |          | #                       | Item #   |
| #                       | Item #   | #                       | Item #   | 6a                      | DMLP900R |
| 6a                      | DMLP605R | 6a                      | DMLP505R |                         |          |
| 6b                      | DMLP505R | 6b                      | DMSP805R |                         |          |

Hide Ray Data

#### RAY DATA

Ray data for Zemax is available for some of the bare LEDs incorporated into these high-powered light sources. This data is provided in a zipped folder that can be downloaded by clicking on the red document

| Item #                | Information File          | Available Ray Files                            | File Size | Click to Download |
|-----------------------|---------------------------|--|-----------|-------------------|
| M365L2                | M365_Info.pdf             | 100,000 Rays and 1 Million Rays                | 27.4 MB   |                   |
| M385L2                | M385_Info.pdf             | 1 Million Rays and 5 Million Rays              | 148 MB    |                   |
| M405L2                | M405_Info.pdf             | 1 Million Rays                                 | 33.1 MB   |                   |
| M450LP1 <sup>a</sup>  | LD_CQAR_20150731_info.pdf | 100,000 Rays, 500,000 Rays, and 5 Million Rays | 123 MB    |                   |
| M455L3 <sup>a,b</sup> | LD_CQ7P_290311_info.pdf   | 100,000 Rays, 500,000 Rays, and 5 Million Rays | 125 MB    |                   |





|        |  |                            |                             |                                  |                      |                                 |                       |            |                     |
|--------|--|----------------------------|-----------------------------|----------------------------------|----------------------|---------------------------------|-----------------------|------------|---------------------|
| 370 nm | (2.5 mW)                                       | -                          | -                           | -                                | -                    | -                               | -                     | -          | -                   |
| 375 nm | -  | M375D2<br>(387 mW<br>Min)  | M375L3<br>(387 mW<br>Min)   | -                                | M375F2<br>(4.23 mW)  | -                               | -                     | -          | -                   |
| 385 nm | -  | M385D1<br>(270 mW<br>Min)  | M385L2<br>(270 mW<br>Min)   | M385L2<br>(90 mW) <sup>d</sup>   | M385F1<br>(10.7 mW)  | -                               | Available<br>(95 mW)  | -          | -                   |
|        |  | M385D2<br>(1650 mW<br>Min) | M385LP1<br>(1650 mW<br>Min) | M385LP1<br>(520 mW) <sup>d</sup> | M385FP1<br>(23.2 mW) |                                 |                       |            |                     |
| 395 nm | -  | M395D3<br>(400 mW<br>Min)  | M395L4<br>(400 mW<br>Min)   | -                                | M395F3<br>(6.8 mW)   | -                               | -                     | -          | -                   |
| 405 nm | LED405E<br>(10 mW)                             | M405D1<br>(410 mW<br>Min)  | M405L2<br>(410 mW<br>Min)   | M405L2<br>(260 mW) <sup>d</sup>  | M405F1<br>(3.7 mW)   | -                               | Available<br>(95 mW)  | DC3100-405 | -                   |
|        |  | M405D2<br>(1500 mW<br>Min) | M405LP1<br>(1500 mW<br>Min) | M405LP1<br>(450 mW) <sup>d</sup> | M405FP1<br>(24.3 mW) |                                 |                       |            |                     |
| 420 nm | -  | M420D2<br>(750 mW<br>Min)  | M420L3<br>(750 mW<br>Min)   | -                                | M420F2<br>(16.2 mW)  | -                               | Available<br>(290 mW) | -          | -                   |
| 450 nm | -  | M450D3<br>(1850 mW<br>Min) | M450LP1<br>(1850 mW<br>Min) | -                                | -                    | -                               | -                     | -          | -                   |
| 455 nm | -  | M455D2<br>(900 mW<br>Min)  | M455L3<br>(900 mW<br>Min)   | M455L3<br>(360 mW) <sup>d</sup>  | M455F1<br>(11.0 mW)  | -                               | Available<br>(310 mW) | -          | -                   |
| 460 nm | -  | -                          | -                           | -                                | -                    | SOLIS-<br>460A(/M)<br>(4175 mW) | -                     | -          | -                   |
| 465 nm | LED465E<br>(20 mW)                             | -                          | -                           | -                                | -                    | -                               | -                     | -          | -                   |
| 470 nm | LED470L<br>(170 mW)                            | M470D2<br>(650 mW<br>Min)  | M470L3<br>(650 mW<br>Min)   | M470L3<br>(250 mW) <sup>d</sup>  | M470F1<br>(10.1 mW)  | -                               | Available<br>(250 mW) | DC3100-470 | LIU470A<br>(253 mW) |
| 490 nm | -  | M490D2<br>(200 mW<br>Min)  | M490L3<br>(200 mW<br>Min)   | -                                | M490F2<br>(2.0 mW)   | -                               | Available<br>(50 mW)  | -          | -                   |
| 505 nm | -  | M505D2<br>(400 mW<br>Min)  | M505L3<br>(400 mW<br>Min)   | M505L3<br>(150 mW) <sup>d</sup>  | M505F1<br>(8.0 mW)   | -                               | Available<br>(170 mW) | -          | -                   |
| 525 nm | LED525E<br>(2.6 mW Max)<br>LED528EHP<br>(7 mW) | -                          | -                           | -                                | -                    | SOLIS-<br>525A(/M)<br>(1650 mW) | -                     | -          | LIU525A<br>(111 mW) |
| 530 nm | -  | M530D2<br>(350 mW<br>Min)  | M530L3<br>(350 mW<br>Min)   | M530L3<br>(130 mW) <sup>d</sup>  | M530F1<br>(5.1 mW)   | -                               | Available<br>(100 mW) | -          | -                   |
| 565 nm | -  | M565D2<br>(880 mW<br>Min)  | M565L3<br>(880 mW<br>Min)   | -                                | M565F1<br>(2.0 mW)   | -                               | Available<br>(106 mW) | -          | -                   |
| 590 nm | LED591E<br>(2 mW)                              | M590D2<br>(160 mW<br>Min)  | M590L3<br>(160 mW<br>Min)   | M590L3<br>(60 mW) <sup>d</sup>   | M590F1<br>(3.2 mW)   | -                               | Available<br>(65 mW)  | -          | LIU590A<br>(109 mW) |
| 595 nm | -  | M595D2<br>(445 mW<br>Min)  | M595L3<br>(445 mW<br>Min)   | -                                | -                    | -                               | -                     | -          | -                   |
| 617 nm | -  | M617D2<br>(600 mW<br>Min)  | M617L3<br>(600 mW<br>Min)   | M617L3<br>(230 mW) <sup>d</sup>  | M617F1<br>(10.8 mW)  | -                               | Available<br>(210 mW) | -          | -                   |
| 623 nm | -  | -                          | -                           | -                                | -                    | SOLIS-<br>623A(/M)<br>(2530 mW) | -                     | -          | -                   |
| 625 nm | -  | M625D2<br>(700 mW<br>Min)  | M625L3<br>(700 mW<br>Min)   | M625L3<br>(270 mW) <sup>d</sup>  | M625F1<br>(10.1 mW)  | -                               | Available<br>(240 mW) | -          | -                   |
| 630 nm | -  | -                          | -                           | -                                | -                    | -                               | -                     | DC3100-630 | LIU630A<br>(208 mW) |
| 635 nm | LED631E<br>(4 mW)                              | -                          | -                           | -                                | -                    | -                               | -                     | -          | -                   |
|        | LED635L<br>(170 mW)                            | -                          | -                           | -                                | -                    | -                               | -                     | -          | -                   |
| 639 nm | LED630E<br>(7.2 mW)                            | -                          | -                           | -                                | -                    | -                               | -                     | -          | -                   |
| 660 nm | -  | M660D2<br>(940 mW<br>Min)  | M660L4<br>(940 mW<br>Min)   | M660L4<br>(400 mW) <sup>d</sup>  | M660F1<br>(14.5 mW)  | -                               | Available<br>(210 mW) | -          | -                   |
| 730 nm | -  | M730D2<br>(515 mW<br>Min)  | M730L4<br>(515 mW<br>Min)   | M730L4<br>(165 mW) <sup>d</sup>  | -                    | -                               | -                     | -          | -                   |

|                                      |  |                           |                           |                                 |                     |                           |   |   |                     |
|--------------------------------------|--|---------------------------|---------------------------|---------------------------------|---------------------|---------------------------|---|---|---------------------|
| 740 nm                               | -  | -                         | -                         | -                               | M740F2<br>(6.0 mW)  | -                         | - | - | -                   |
| 780 nm                               | LED780E<br>(18 mW)                                       | M780D2<br>(200 mW<br>Min) | M780L3<br>(200 mW<br>Min) | M780L3<br>(130 mW) <sup>d</sup> | M780F2<br>(7.5 mW)  | -                         | - | - | LIU780A<br>(315 mW) |
| 810 nm                               | -  | M810D2<br>(325 mW<br>Min) | M810L3<br>(325 mW<br>Min) | M810L3<br>(210 mW) <sup>d</sup> | M810F2<br>(6.5 mW)  | -                         | - | - | -                   |
| 850 nm                               | LED851W<br>(8 mW)<br>LED851L<br>(13 mW)                  | M850D2<br>(900 mW<br>Min) | M850L3<br>(900 mW<br>Min) | M850L3<br>(330 mW) <sup>d</sup> | M850F2<br>(13.4 mW) | -                         | - | - | LIU850A<br>(322 mW) |
| 870 nm                               | LED870E<br>(22 mW)                                       | -                         | -                         | -                               | -                   | -                         | - | - | -                   |
| 880 nm                               | -  | M880D2<br>(300 mW<br>Min) | M880L3<br>(300 mW<br>Min) | -                               | M880F2<br>(3.4 mW)  | -                         | - | - | -                   |
| 910 nm                               | LED910E<br>(12 mW)                                       | -                         | -                         | -                               | -                   | -                         | - | - | -                   |
| 940 nm                               | LED940E<br>(18 mW)                                       | M940D2<br>(800 mW<br>Min) | M940L3<br>(800 mW<br>Min) | M940L3<br>(320 mW) <sup>d</sup> | M940F1<br>(6.5 mW)  | -                         | - | - | -                   |
| 970 nm                               | -  | M970D2<br>(35 mW Min)     | M970L3<br>(35 mW Min)     | -                               | M970F2<br>(0.3 mW)  | -                         | - | - | -                   |
| 1050 nm                              | LED1050E<br>(2.5 mW)                                     | M1050D1<br>(50 mW Min)    | M1050L2<br>(50 mW Min)    | -                               | M1050F1<br>(1.4 mW) | -                         | - | - | -                   |
| 1070 nm                              | LED1070E<br>(7.5 mW)                                     | -                         | -                         | -                               | -                   | -                         | - | - | -                   |
| 1200 nm                              | LED1200E<br>(2.5 mW)                                     | M1200D2<br>(30 mW Min)    | M1200L3<br>(30 mW Min)    | -                               | -                   | -                         | - | - | -                   |
| 1300 nm                              | LED1300E<br>(2 mW)                                       | M1300D2<br>(25 mW Min)    | M1300L3<br>(25 mW Min)    | -                               | -                   | -                         | - | - | -                   |
| 1450 nm                              | LED1450E<br>(2 mW)                                       | M1450D2<br>(31 mW Min)    | M1450L3<br>(31 mW Min)    | -                               | -                   | -                         | - | - | -                   |
| 1550 nm                              | LED1550E<br>(2 mW)                                       | M1550D2<br>(31 mW Min)    | M1550L3<br>(31 mW Min)    | -                               | -                   | -                         | - | - | -                   |
| 1650 nm                              | LED1600P<br>(1.2 mW)                                     | -                         | -                         | -                               | -                   | -                         | - | - | -                   |
| 1750 nm                              | LED1700P<br>(1.2 mW<br>Quasi-CW,<br>30 mW<br>Pulsed)     | -                         | -                         | -                               | -                   | -                         | - | - | -                   |
| 1850 nm                              | LED1800P<br>(0.9 mW<br>Quasi-CW,<br>20 mW<br>Pulsed)     | -                         | -                         | -                               | -                   | -                         | - | - | -                   |
| 1950 nm                              | LED1900P<br>(1.0 mW<br>Quasi-CW,<br>25 mW<br>Pulsed)     | -                         | -                         | -                               | -                   | -                         | - | - | -                   |
| 2050 nm                              | LED2050P<br>(1.1 mW<br>Quasi-CW,<br>28 mW<br>Pulsed)     | -                         | -                         | -                               | -                   | -                         | - | - | -                   |
| 2350 nm                              | LED2350P<br>(0.8 mW<br>Quasi-CW,<br>16 mW<br>Pulsed)     | -                         | -                         | -                               | -                   | -                         | - | - | -                   |
| 4200 nm                              | LED4300P<br>(0.01 mW<br>Quasi-CW,<br>0.2 mW<br>Pulsed)   | -                         | -                         | -                               | -                   | -                         | - | - | -                   |
| 4500 nm                              | LED4600P<br>(0.006 mW<br>Quasi-CW,<br>0.12 mW<br>Pulsed) | -                         | -                         | -                               | -                   | -                         | - | - | -                   |
| 467.5 nm,<br>525 nm,<br>and 627.5 nm | LEDRGBE<br>(5.8 mW,<br>6.2 mW,<br>and 3.1 mW)            | -                         | -                         | -                               | -                   | -                         | - | - | -                   |
| 470 - 850 nm                         | -  | MBB1D1<br>(70 mW Min)     | MBB1L3<br>(70 mW Min)     | -                               | MBB1F1<br>(1.2 mW)  | -                         | - | - | -                   |
| 6500 K<br>(Cold White)               | LEDWE-15<br>(13 mW)                                      | MCWHD2<br>(800 mW<br>Min) | MCWHL5<br>(800 mW<br>Min) | MCWHL5<br>(320 mW) <sup>d</sup> | -                   | SOLIS-1A(/M)<br>(3070 mW) | - | - | LIUCWHA<br>(250 mW) |
| 5600 K                               |  | -                         | -                         | -                               | MCWHF1              | -                         | - | - | -                   |

|                        |   |                           |                           |   |                    |  |   |   |   |
|------------------------|---|---------------------------|---------------------------|---|--------------------|--|---|---|---|
| (Cold White)           |   | -                         | -                         | - | (7.0 mW)           |  | - | - |   |
| 3000 K<br>(Warm White) | - | MWWHD1<br>(500 mW<br>Min) | MWWHL3<br>(500 mW<br>Min) | - | MWWHF1<br>(7.0 mW) |  | - | - | - |

- These Collimated LEDs are compatible with the standard and epi-illumination ports on the following microscopes: Olympus BX/IX (Item # Suffix: -C1), Leica DMI (Item # Suffix: -C2), Zeiss Axioskop (Item # Suffix: -C4), and Nikon Eclipse (Bayonet Mount, Item # Suffix: -C5).
- Typical power when used with MM Fiber with Ø400 µm core, 0.39 NA.
- Our LED4D 4-Wavelength LED Source is available with select combinations of the LEDs at these wavelengths.
- Typical power for LEDs with the Leica DMI collimation package (Item # Suffix: -C2).

[Hide Mounted LEDs with EEPROM \(>1000 mW Output Power\)](#)

**Mounted LEDs with EEPROM (>1000 mW Output Power)**



- ▶ Mounted LEDs with >1000 mW (Min) Output
- ▶ Integrated EEPROM for Automated LED Settings
- ▶ >10,000 Hours Lifetime (Except M450LP1; See *Specs* and *Stability* Tab for Details)
- ▶ Integrated Large Heat Sink for Optimized Thermal Management
- ▶ Output can be Modulated with Suitable Controller (See the *LED Drivers* Tab)
- ▶ Compatible with 30 mm Cage System and SM1 Lens Tubes
- ▶ Cable Length: 2 m



Click to Enlarge  
M385LP1 LED Inserted  
into CP02 Cage Plate  
and Mounted Using Ø6  
mm Cage Rods



Click to Enlarge  
M385LP1 LED with  
SM1L03 Lens Tube  
Post Mounted Using an  
SM1RC/M Slip Ring

These LEDs with output powers greater than 1000 mW are mounted to the end of a larger heat sink in a Ø57.0 mm plastic housing. The larger heat sink is necessary in order to dissipate the increased heat output of these diodes. The heat sink housing comes equipped with internal SM1 (1.035"-40) threads and four 4-40 tapped holes for compatibility with Thorlabs' SM1 lens tubes and 30 mm cage systems, respectively.

Please note that these mounted LEDs with output wavelengths of 365 nm, 385 nm, or 405 nm radiate intense UV light during operation. Precautions must be taken to prevent looking directly at the UV light and UV light protective glasses must be worn to avoid eye damage. Exposure of the skin and other body parts to the UV light should be avoided. These LEDs are not intended for use in household illumination applications.

| Part Number | Description  | Price    | Availability |
|-------------|--|----------|--------------|
| M365LP1     | NEW! UV (365 nm) Mounted LED, 1400 mA, 1150 mW (Min)         | \$416.67 | Today        |
| M385LP1     | NEW! UV (385 nm) Mounted LED, 1400 mA, 1650 mW (Min)         | \$416.67 | Today        |
| M405LP1     | NEW! UV (405 nm) Mounted LED, 1400 mA, 1500 mW (Min)         | \$416.67 | Today        |
| M450LP1     | NEW! Royal Blue (450 nm) Mounted LED, 2000 mA, 1850 mW (Min) | \$294.44 | Today        |

[Hide Mounted LEDs with EEPROM \(<1000 mW Output Power\)](#)

**Mounted LEDs with EEPROM (<1000 mW Output Power)**



- ▶ Integrated EEPROM for Automated LED Settings
- ▶ Long Lifetimes (See *Specs* and *Stability* Tab for Details)
  - >10,000 Hours for LEDs with a Nominal Wavelength of ≥365 nm
  - >500 Hour Lifetime for LEDs with a Nominal Wavelength of ≤340 nm
- ▶ Stable Output Intensity by Optimized Thermal Management
- ▶ Output can be Modulated with Suitable Controller (See the *LED Drivers* Tab)
- ▶ Compatible with Thorlabs' SM1 Lens Tubes
- ▶ Fits Inside a 30 mm Cage System
- ▶ Cable Length: 2 m



Click to Enlarge  
MWWHL3 LED  
Mounted in  
an SM1RC Slip Ring

These LEDs with output powers less than 1000 mW are mounted to the end of a Ø30.5 mm heatsink for heat dissipation and thermal stability. The heat sink is equipped with internal SM1 (1.035"-40) threads for compatibility with Thorlabs' SM1 lens tubes.

Please note that our LEDs with wavelengths from 280 nm to 420 nm radiate intense UV light during operation. Precautions must be taken to prevent looking directly at the UV light and UV light protective glasses must be worn to avoid eye damage. Exposure of the skin and other body parts to the UV light should be avoided. Mounted LEDs are not intended for use in household illumination applications.


| Part Number | Description   | Price      | Availability |
|-------------|---|------------|--------------|
| M265L3      | NEW! Deep UV (265 nm) Mounted LED, 350 mA, 10 mW (Min)          | \$1,230.00 | Today        |
| M280L3      | Deep UV (280 nm) Mounted LED, 350 mA, 25 mW (Min)               | \$1,230.00 | Today        |
| M310L3      | Deep UV (310 nm) Mounted LED, 350 mA, 25 mW (Min)               | \$1,230.00 | Today        |
| M340L3      | Customer Inspired!UV (340 nm) Mounted LED, 80 mA, 10 mW (Min)   | \$1,100.00 | Today        |
| M365L2      | UV (365 nm) Mounted LED, 700 mA, 190 mW (Min)                   | \$469.00   | Today        |
| M375L3      | Customer Inspired!UV (375 nm) Mounted LED, 700 mA, 387 mW (Min) | \$247.00   | Today        |
| M385L2      | UV (385 nm) Mounted LED, 700 mA, 270 mW (Min)                   | \$469.00   | Today        |
| M395L4      | UV (395 nm) Mounted LED, 500 mA, 400 mW (Min)                   | \$283.00   | Today        |
| M405L2      | UV (405 nm) Mounted LED, 1000 mA, 410 mW (Min)                  | \$469.00   | Today        |
| M420L3      | Violet (420 nm) Mounted LED, 1000 mA, 750 mW, (Min)             | \$283.00   | Today        |
| M455L3      | Royal Blue (455 nm) Mounted LED, 1000 mA, 900 mW (Min)          | \$268.00   | Today        |
| M470L3      | Blue (470 nm) Mounted LED, 1000 mA, 650 mW (Min)                | \$268.00   | Today        |
| M490L3      | Blue (490 nm) Mounted LED, 350 mA, 200 mW (Min)                 | \$268.00   | Today        |
| M505L3      | Cyan (505 nm) Mounted LED, 1000 mA, 400 mW (Min)                | \$268.00   | Today        |

|         |   |          |           |
|---------|---|----------|-----------|
| M530L3  | Green (530 nm) Mounted LED, 1000 mA, 350 mW (Min)                   | \$268.00 | Today     |
| M565L3  | Lime (565 nm) Mounted LED, 1000 mA, 880 mW (Min)                    | \$211.00 | Today     |
| M590L3  | Amber (590 nm) Mounted LED, 1000 mA, 160 mW (Min)                   | \$193.00 | Today     |
| M595L3  | Phosphor-Converted Amber (595 nm) Mounted LED, 700 mA, 445 mW (Min) | \$211.00 | Today     |
| M617L3  | Orange (617 nm) Mounted LED, 1000 mA, 600 mW (Min)                  | \$193.00 | Today     |
| M625L3  | Red (625 nm) Mounted LED, 1000 mA, 700 mW (Min)                     | \$193.00 | Today     |
| M660L4  | NEW! Deep Red (660 nm) Mounted LED, 1200 mA, 940 mW (Min)           | \$211.11 | Lead Time |
| M730L4  | Far Red (730 nm) Mounted LED, 1000 mA, 515 mW (Min)                 | \$190.00 | Today     |
| M780L3  | IR (780 nm) Mounted LED, 800 mA, 200 mW (Min)                       | \$211.00 | Today     |
| M810L3  | IR (810 nm) Mounted LED, 500 mA, 325 mW (Min)                       | \$196.67 | Today     |
| M850L3  | IR (850 nm) Mounted LED, 1000 mA, 900 mW (Min)                      | \$211.00 | Today     |
| M880L3  | IR (880 nm) Mounted LED, 1000 mA, 300 mW (Min)                      | \$211.00 | Today     |
| M940L3  | IR (940 nm) Mounted LED, 1000 mA, 800 mW (Min)                      | \$211.00 | Today     |
| M970L3  | IR (970 nm) Mounted LED, 600 mA, 35 mW (Min)                        | \$211.00 | Today     |
| M1050L2 | Customer Inspired!!IR (1050 nm) Mounted LED, 700 mA, 50 mW (Min)    | \$227.00 | Today     |
| M1200L3 | Customer Inspired!!IR (1200 nm) Mounted LED, 700 mA, 30 mW (Min)    | \$282.00 | Today     |
| M1300L3 | Customer Inspired!!IR (1300 nm) Mounted LED, 500 mA, 25 mW (Min)    | \$282.00 | Today     |
| M1450L3 | IR (1450 nm) Mounted LED, 700 mA, 31 mW (Min)                       | \$282.00 | Today     |
| M1550L3 | Customer Inspired!!IR (1550 nm) Mounted LED, 700 mA, 31 mW (Min)    | \$282.00 | Today     |
| MBB1L3  | Broadband (470 - 850 nm) Mounted LED, 500 mA, 70 mW (Min)           | \$494.00 | Today     |
| MWWHL3  | Warm White Mounted LED, 1000 mA, 500 mW (Min)                       | \$193.00 | Today     |
| MCWHL5  | Cold White Mounted LED, 1000 mA, 800 mW (Min)                       | \$193.00 | Today     |

[Hide Collimation Adapters with AR-Coated Ø50 mm Aspheric Lenses](#)

### Collimation Adapters with AR-Coated Ø50 mm Aspheric Lenses

- ▶ AR-Coated Aspheric Lens with Low f/# (Approximately 0.8)
- ▶ Compatible with Select Leica, Nikon, Olympus, or Zeiss Microscopes
- ▶ Easily Adjust Beam Collimation / Focus
- ▶ Requires SM2T2 Coupler and SM1A2 Adapter (Each Sold Separately) when Used with the LEDs Above

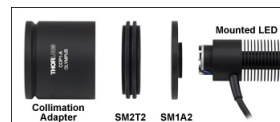

[Click for Details](#)

Thorlabs offers collimation adapters with Ø2" AR-coated aspheric condenser lenses (EFL: 40 mm) for collimating the output from the mounted LEDs sold above. Two AR coating ranges (350 - 700 nm and 650 - 1050 nm) and four different collimator housings are available. Each housing is designed to mate to the illumination port on selected Olympus\*, Leica, Nikon, or Zeiss microscopes. Compatible microscopes are listed in the Collimation Adapter Selection Guide table below.

Installation of Collimation Adapter to Mounted LED Using SM2T2 and SM1A2

Using an adapter with a substrate or AR coating that does not transmit the wavelength of your LED is not recommended. The deep UV LEDs (M280L3, M310L3, M340L3) require a lens fabricated from UV Fused Silica, since many standard varieties of glass do not transmit below 350 nm. IR LEDs that emit beyond 1050 nm (M1200L3, M1300L3, M1450L3, M1550L3) can be collimated using an uncoated condenser lens; the ACL5040U is an uncoated version of the Ø2" lenses used in the collimation packages below that has a wavelength range of 380 - 2100 nm. Alternatively, the ACP collimation adapter below accepts user-supplied Ø1" collimation optics and includes a thread adapter that converts the internally M34 x 0.5 threaded output to our SM2 (2.035" -40) thread standard. See the *Collimation* tab above for more information on collimation options.

The LED sources described above can be fitted to the collimators by using an SM2T2 Coupler and SM1A2 Adapter (not included) as shown in the image at right. This assembly can be easily adapted to different LED sources by unscrewing the LED housing.

If compatibility with SM1 (1.035"-40) threading is preferable to compatibility with a microscope for your application, the ACP collimation adapters available below can be used with the SM1A38 M34 x 0.5 to SM1 thread adapter.

\* Please note that due to the optical design of the transmitted lamphouse port of the BX and IX microscopes it may be necessary to purchase a separate adapter which is available from Olympus.

| Collimation Adapter Selection Guide |                   |             |   |   |  |   |
|-------------------------------------|-------------------|-------------|---|---|--|---|
| Compatible Microscopes              |                   |             | Olympus BX & IX <sup>a</sup>  | Leica DMI   | Zeiss Axioskop   | Nikon Eclipse   |
| AR Coating Range of Condenser Lens  | Lens Focal Length | Lens Item # |  |  |  |  |
|                                     |                   |             | Click to Enlarge  | Click to Enlarge  | Click to Enlarge   | Click to Enlarge  |
|                                     |                   |             | 350 - 700 nm  | 40.0 mm   | ACL5040U-A   | COP1-A  |
| 650 - 1050 nm                       | 40.0 mm           | ACL5040U-B  | COP1-B  | COP2-B  | COP4-B   | COP5-B  |

- Please note that due to the optical design of the transmitted lamphouse port of the BX and IX microscopes it may be necessary to purchase a separate adapter which is available from Olympus.

| Part Number | Description  | Price    | Availability |
|-------------|--|----------|--------------|
| COP1-A      | Collimation Adapter for Olympus BX & IX, AR Coating: 350 - 700 nm  | \$181.00 | Today        |
| COP1-B      | Collimation Adapter for Olympus BX & IX, AR Coating: 650 - 1050 nm | \$211.00 | Today        |
| COP2-A      | Collimation Adapter for Leica DMI, AR Coating: 350 - 700 nm        | \$181.00 | Today        |
| COP2-B      | Collimation Adapter for Leica DMI, AR Coating: 650 - 1050 nm       | \$211.00 | Today        |
| COP4-A      | Collimation Adapter for Zeiss Axioskop, AR Coating: 350 - 700 nm   | \$181.00 | Today        |
| COP4-B      | Collimation Adapter for Zeiss Axioskop, AR Coating: 650 - 1050 nm  | \$211.00 | Today        |
| COP5-A      | Collimation Adapter for Nikon Eclipse, AR Coating: 350 - 700 nm    | \$214.00 | Today        |
| COP5-B      | Collimation Adapter for Nikon Eclipse, AR Coating: 650 - 1050 nm   | \$249.00 | Today        |

|       |  |         |       |
|-------|--|---------|-------|
| SM1A2 | Adapter with External SM1 Threads and Internal SM2 Threads | \$24.00 | Today |
| SM2T2 | SM2 (2.035"-40) Coupler, External Threads                  | \$34.00 | Today |

[Hide Collimation Adapters for Ø1" or Ø25 mm Optics](#)

### Collimation Adapters for Ø1" or Ø25 mm Optics



- ▶ Integrate a Ø1" or Ø25 mm Collimation Optic with Thorlabs' High-Power Mounted LEDs
- ▶ Three Versions Available:
  - ▶ ACP2520-A Includes Ø25 mm AR-Coated Lens for 350 - 700 nm
  - ▶ ACP2520-B Includes Ø25 mm AR-Coated Lens for 650 - 1050 nm
  - ▶ ACP Mechanical Housing without Lens Accepts User-Supplied Ø1" or Ø25 mm Optics
- ▶ Adjust and Set Lens Position via Rotating Ring with Locking Setscrew
- ▶ Removable External M34 x 0.5 to Internal SM2 (2.035"-40) Thread Adapter
- ▶ Compatible with Thorlabs' Internally SM2-Threaded Microscope Port Adapters

| Item #           | Included Lens | AR Coating Wavelength Range | Lens Focal Length |
|------------------|---------------|-----------------------------|-------------------|
| ACP2520-A        | ACL2520U-A    | 350 - 700 nm                | 20.0 mm           |
| ACP2520-B        | ACL2520U-B    | 650 - 1050 nm               | 20.0 mm           |
| ACP <sup>a</sup> | -             | -                           | -                 |

- The ACP does not include a collimation optic, allowing user-supplied optics to be integrated with Thorlabs' high-power mounted LEDs.



Click to Enlarge  
The ACP2520-A  
Installed on a High-  
Power LED

These adapters allow a Ø1" or Ø25 mm collimation optic to be integrated with the mounted LEDs sold above. The input of each housing is externally SM1 (1.035"-40) threaded to mate with the housing of our high-power LEDs, while the output has internal M34 x 0.5 threads. An external M34 x 0.5 to internal SM2 adapter is included so that these collimation adapters can be used with our SM2-threaded microscope adapters or other SM2-threaded components. Thorlabs also offers an M34 x 0.5 external to SM1 (1.035"-40) thread adapter (Item # SM1A38) that can be used to integrate the collimation adapters with SM1-threaded components.

The ACP2520-A and ACP2520-B include a removable AR-coated aspheric condenser lenses for 350 - 700 nm and 650 - 1050 nm, respectively. The ACP does not include a collimation optic, allowing user-supplied optics to be integrated with Thorlabs' high-power mounted LEDs. The lens is mounted in an inner carriage that can be translated up to 11 mm (0.43") along the Z-axis by turning the adjustment ring (engraved with the item # in the photos to the left) and locked into position by turning the locking screw on the side of the adjustment ring with a 2 mm (5/64") hex key. Lines with a 2 mm spacing are engraved on the housing as a rough guide for how far the carriage has been translated. Each collimation adapter also includes an extra-thick SM1-threaded retaining ring that allows an aspheric condenser lens to be mounted with our SPW602 spanner wrench.

Each collimation adapter also includes an extra-thick SM1-threaded retaining ring that allows an aspheric condenser lens to be mounted with our SPW602 spanner wrench.

#### Inserting or Removing Optics

To insert or remove an optic in these collimation adapters, use the adjustment ring to translate the inner carriage to the M34 x 0.5-threaded end of the housing. Remove the included SM1-threaded retaining ring using a spanner wrench. If there is a lens installed already, remove it from the carriage. Insert another Ø1" or Ø25 mm optic into the carriage, and use the retaining ring to secure it.

Using an adapter with a substrate or AR coating that does not transmit the wavelength of your LED is not recommended. The deep UV LEDs (M280L3, M310L3, M340L3) require a lens fabricated from UV Fused Silica, since many standard varieties of glass do not transmit below 350 nm. IR LEDs that emit beyond 1050 nm (M1200L3, M1300L3, M1450L3, M1550L3) can be collimated using an uncoated condenser lens, such as the Ø25 mm ACL2520U, that has a wavelength range of 380 - 2100 nm.

| Part Number | Description  | Price    | Availability |
|-------------|--|----------|--------------|
| ACP2520-A   | Adjustable Collimation Adapter with AR-Coated Lens for 350 - 700 nm  | \$194.44 | 3-5 Days     |
| ACP2520-B   | Adjustable Collimation Adapter with AR-Coated Lens for 650 - 1050 nm | \$194.44 | Today        |
| ACP         | Adjustable Collimation Adapter for Ø1" or Ø25 mm Optic               | \$175.56 | Today        |

[Hide Mounted LED Mating Connector](#)

### Mounted LED Mating Connector



- ▶ Pico (M8) Receptacle
- ▶ Female 4-Pin for Front Mounting
- ▶ 0.5 m Long, 24 AWG Wires
- ▶ M8 x 0.5 Panel Mount Thread
- ▶ IP 67 and NEMA 6P Rated

The CON8ML-4 connector can be used to mate mounted LEDs featured on this page to user-supplied power supplies. We also offer a male 4-Pin M8 connector cable (Item # CAB-LEDD1).

| Pin | Color | Specification |
|-----|-------|---------------|
| 1   | Brown | LED Anode     |
| 2   | White | LED Cathode   |
| 3   | Black | EEPROM GND    |
| 4   | Blue  | EEPROM IO     |



| Part Number | Description                                    | Price   | Availability |
|-------------|--|---------|--------------|
| CON8ML-4    | 4-Pin Female Mating Connector for Mounted LEDs | \$30.00 | Today        |