

CHAPTERS

Tables/
Breadboards

Mechanics

Optomechanic
Devices

Kits

Lab Supplies

SECTIONS

Mini Series

Ø1/2" Post
Assemblies

Ø1" Post
Assemblies

Ø1.5" Post
Assemblies

Mounting/Angle
Brackets

Lab Platforms

Instrument Shelves

Lens Tubes

Cage Systems

Optical Rails

Mirror Mounts

Kinematic Mounts

Fixed Mounts

Translation Mounts

Rotation Mounts

Filter Mounts

V-Mounts

Iris/Apertures

Adapters

SM1 Series of Stackable Lens Tubes for Ø1" Optics



SM1L20



SM1L03



SM1L05




SM1L10




SM1L30

These stackable Ø1.00" lens tubes provide a fast and easy means for building complex optical assemblies. Multiple optics can be inserted into a single tube using multiple retaining rings (SM1RR); simply insert the first ring at the depth required, add your optic, and thread a second retaining ring to secure the optic.

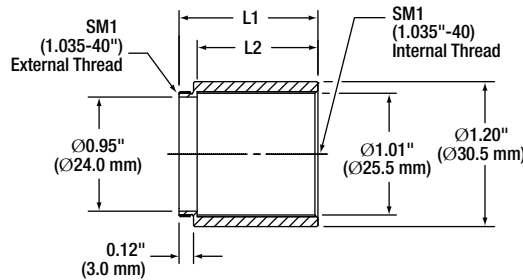

Also Available in
Ø1/2", Ø30 mm,
Ø2.00", and Ø3.00"

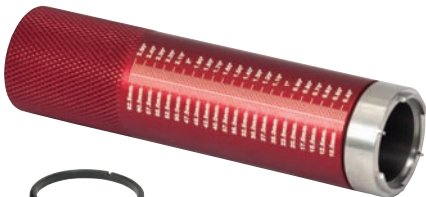

SPW602
Spanner Wrench

ITEM #	\$	£	€	RMB	L1	L2	DESCRIPTION
SM1L03*	\$ 12.16	£ 8.76	€ 10.58	¥ 96.92	0.45"	0.30"	SM1 Lens Tube, 0.3" Deep
SM1L05*	\$ 12.59	£ 9.06	€ 10.95	¥ 100.34	0.65"	0.50"	SM1 Lens Tube, 1/2" Deep
SM1L10*	\$ 14.25	£ 10.26	€ 12.40	¥ 113.57	1.15"	1.00"	SM1 Lens Tube, 1.0" Deep
SM1L20*	\$ 16.50	£ 11.88	€ 14.36	¥ 131.51	2.15"	2.00"	SM1 Lens Tube, 2.0" Deep
SM1L30*	\$ 25.75	£ 18.54	€ 22.40	¥ 205.23	3.15"	3.00"	SM1 Lens Tube, 3.0" Deep
SM1RR	\$ 4.50	£ 3.24	€ 3.92	¥ 35.87	SM1 (1.035"-40) Threaded Retaining Ring		
SPW602	\$ 26.00	£ 18.72	€ 22.62	¥ 207.22	Spanner Wrench for SM1RR		

*One SM1RR Retaining Ring Included

SM1 Series Lens Tubes Without External Threads



SM1M05




SM1M10



SM1M20




SM1RR

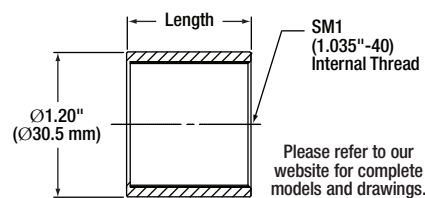

Also Available in
Ø1/2" and
Ø2.00"

Please refer to our
website for complete
models and drawings.

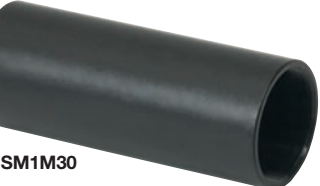

SM1M30

Unlike our Ø1.00" stackable lens tubes featured above, these internally threaded lens tubes do not feature a lip at one end nor do they have external SM1 threading. Consequently, optics housed inside can be accessed from both ends. To join these lens tubes end-to-end, use an externally-threaded lens tube coupler (see page 137) or a lens tube without internal threads (see following page). Two SM1RR retaining rings are included with each lens tube.

ITEM #	\$	£	€	RMB	DESCRIPTION
SM1M05*	\$ 12.00	£ 8.64	€ 10.44	¥ 95.64	SM1 Lens Tube, Without External Threads, 1/2" Long
SM1M10*	\$ 13.25	£ 9.54	€ 11.53	¥ 105.60	SM1 Lens Tube, Without External Threads, 1.00" Long
SM1M20*	\$ 15.50	£ 11.16	€ 13.49	¥ 123.54	SM1 Lens Tube, Without External Threads, 2.00" Long
SM1M30*	\$ 16.75	£ 12.06	€ 14.57	¥ 133.50	SM1 Lens Tube, Without External Threads, 3.00" Long
SM1RR	\$ 4.50	£ 3.24	€ 3.92	¥ 35.87	SM1 (1.035"-40) Threaded Retaining Ring

*Two SM1RR Retaining Rings Included