

Ø1.0" Pillar Post C-Clamps

The RM1A and RM1B C-Clamps for Ø1" posts have two mounting points located 90° apart for perpendicular mounting of components. One of the mounting points has an #8-32 (M4 x 0.7) tapped mounting hole while the other contains an

#8 (M4) counterbored through hole. The C-clamp design allows the clamp to be slipped on and off the post easily by simply loosening or tightening the brass-tipped (non-marring) setscrew that locks the C-clamp to the Ø1" post.

- Slip-On C-Clamp
- Two Mounting Holes: One Counterbore and One Tap



ITEM#	METRIC ITEM#	\$	£	€	RMB	DESCRIPTION
RM1A	RM1A/M	\$ 25.20	£ 17.50	€ 22,40	¥212.80	#8-32 (M4) Pillar Post C-Clamp
RM1B	RM1B/M	\$ 25.20	£ 17.50	€ 22,40	¥212.80	1/4"-20 (M6) Pillar Post C-Clamp

Ø1.0" Pillar Post Angle Clamp

The RM1C Post Angle Clamp has one mounting point with a rotating, threaded mounting stud allowing an affixed optomechanical component or structure to be rotated while maintaining solid contact with the rigid mounting post. The C-clamp design allows the clamp to be slipped on and off the post easily. Simply loosen or tighten the brass-tipped (non-marring) setscrew that locks the C-clamp to the Ø1" post.




ITEM#	METRIC ITEM#	\$	£	€	RMB	DESCRIPTION
RM1C	RM1C/M	\$ 47.40	£ 32.90	€ 42,10	¥ 400.30	Pillar Post Angle Clamp

Ø1.0" Pillar Post, Adjustable-Height C-Clamp

The RM1D and RM1E Variable Offset C-Clamps for Ø1.0" posts have a variable-position mounting surface allowing for adjustability over a range of 0.4" (10 mm) as well as a fixed position tapped mounting hole on the back side of the clamp. The variable-position mounting surface has either a tapped or counterbored

through hole (model dependent) and can be locked into place using a thumbscrew. The C-clamp design allows the clamp to be slipped on and off the post easily. Simply loosen or tighten the brass-tipped (non-marring) setscrew that locks the C-clamp to the Ø1.0" post.

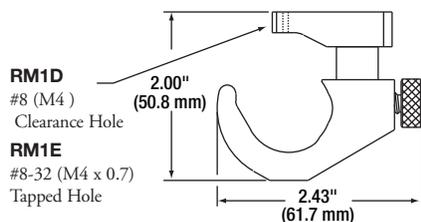

RM1D
#8 (M4)
Clearance Hole

RM1E
#8-32 (M4 x 0.7)
Tapped Hole

RM1D
Clearance Hole for #8 (M4) Screws

RM1E
Tapped Mounting Hole for #8-32 (M4 x 0.7) Screws

Allows 0.40" (10.2 mm) of Vertical Adjustment



Please refer to our website for complete models and drawings.

ITEM#	METRIC ITEM#	\$	£	€	RMB	DESCRIPTION
RM1D	RM1D/M	\$ 50.40	£ 35.00	€ 44,80	¥425.60	Adjustable-Height C-Clamp, #8 (M4) Clearance
RM1E	RM1E/M	\$ 50.40	£ 35.00	€ 44,80	¥425.60	Adjustable-Height C-Clamp, Tapped #8-32 (M4 x 0.7)